Features Denta-Pro

1 Electric micromotor with LED light (0-40'000

Internal spray

Clockwise and anticlockwise rotation

1 Electric scaler with LED light

Automatic instrument recognition

Continuously regulation of the motor and scaler with rotary knob

Continuously regulation of the water quantity with rotary knob

1 Spray water tank: 1 litre

Water filte

1 Foot control

Each instrument is individually adjustable

Features Minimax-Vac

1 High-performance suction hose 300 I/min

1 Rotary knob for continuously variable regulation

Unique BPR Swiss technology

With adapter for saliva ejector tips and oral-drycup

1 Suction fluid container: 1 litre (automatic overflow protection)

1 DENTA-BAG carrying bag

Technical Data

Weight: 2.3 kg (weight including DENTA-BAG and foot pedal: 3.6 kg)

Dimensions DENTA-PRO: 275mm × 130mm × 130mm

Weight: 3.9 kg (weight including DENTA-BAG 4.5 kg)

Dimensions MINIMAX-VAC: 190 mm × 140 mm × 285 mm

Dimensions DENTA-BAG carrying bag: 300mm × 275mm × 275mm

Warranty: 2 years

Options: Contra-angles, MINIMAX-VAC, DENTA-TROLLEY, DENTA-CASE etc.

Denta-Pro & Minimax-Vac


Small and light power couple

The DENTA-PRO is your reliable, efficient and user-friendly partner in mobile dentistry.

To activate the DENTA-PRO simply insert the power plug into an outlet. Despite the very low weight of 2.3 kg the instruments of the DENTA-PRO have complete performance and can be used in continuous operation. In addition to the reliability and performance, the automatic instrument recognition and the internal spray additionally increase the efficiency of the DENTA-PRO.

The MINIMAX-VAC sets new standards in portable plug-and-play suction units.

To start the MINIMAX-VAC it only needs to be plugged into mains power. With a weight of only 3.9 kilograms (!), the MINIMAX-VAC, with over 300I/min (!), achieves a comparable suction to a permanently installed treatment unit in a dental office. The suction can be seamlessly regulateted any time. Due to this pioneering key data MINIMAX-VAC can be used in many areas such as mobile dentistry, in the dental office, or in numerous other areas of medicine.

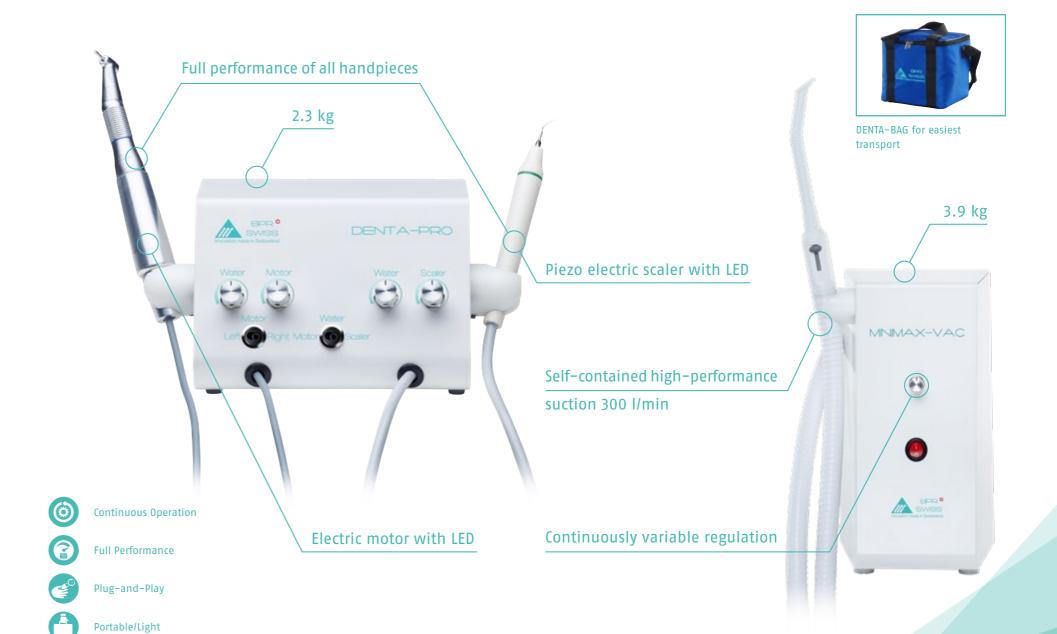

Plug-and-play portable dental units

Denta-Pro & Minimax-Vac