Features

Autonomously functioning cart with integrated oil-free compressor, silent suction system and fresh water tanks

High quality stainless steel housing

1 Three-way syringe (air/water/spray) with sterilizable tips

1 Saliva / surgical suction hose

1 High-performance suction hose 300 I/mi

Continuously adjustable

Unique BPR Swiss technology

2 Fresh water tanks: 2 × 2 litres

2 Suction fluid containers: 2 × 2 litres (automatic overflow protection)

Integrated suction with amalgam separator

Integrated oil-free compressor: Max. operating pressure of 4 bar

Technical Data

Operating pressure: 3.0-4.0 bar

Total weight: 55 kg

Dimensions: 355mm × 690mm × 895mr


Warranty: 2 years

Options: Curing light

Pro-Vac Uniquely convenient flexibility

The mobile and fully stand-alone suction cart PRO-VAC supports you with unlimited performance, cost savings and full flexibility in everyday dental practice and clinic planning.

The PRO-VAC makes your work easier with a three-way syringe and an impressive suction capacity (300 l/min) that outperforms many fixed units. Despite outstanding performance, the oil-free compressor and integrated suction are very quiet and reliable. Dental clinic workflows are simplified by the already integrated amalgam separator. The compact design of the PRO-VAC helps the dental team to efficiently maintain optimal hygiene. Maximum flexibility combined with highest performance and comfort.

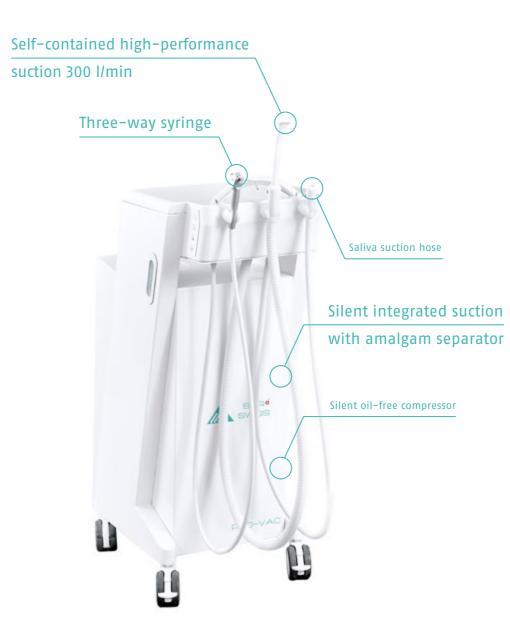



Full Performance

Plug-and-Play

Mobile/Rollable