Features

1 Three-way syringe (air/water/spray) with sterilisable tips (hose and tips detachable)

1 Electric micromotor with LED light (0-40,000)

1 Electric scaler with LED light

Instrument hoses detachable

Spray adjustment regulators for each instrument

Integrated spray water system for syringe and instruments (check valve prevents the back flow of water)

1 Spray water tank: 0.75 litres

1 High-performance suction hose 300 l/min. Continuously adjustable. Unique BPR Swiss technology

Adapter for saliva ejector tips and oral-dry-cup

1 Suction fluid container: 1 litre (automatic overflow protection

Pressure regulator-filter

Sturdy, air- and water-tight case, will not scratch, dent or corrode

Detachable front cover

Disc type foot control

Integrated silent compressor – oil-free: Max. Pressure: 4 ba

DENTA-TROLLEY with drawers for instruments and consumables

Technical Data

Operating pressure: 3-4 bar

Each instrument is individually adjustable

Weight SMART-PORT: 12 kg

Weight DENIA-IRULLEY: 8 kg (space for instruments and consumable

Dimensions SMART-PORT: 380mm × 205mm × 530mm

Dimensions DENTA-TROLLEY: 410mm × 400mm × 550mm

Warranty: 2 years

Options: Contra-angles, air polisher module, additional electrical motor, DENTA-CASE etc.

Smart-Port


The lightness of being

The cleverly designed SMART-PORT integrates the most important functions and operations of a dental treatment unit in a single 12-kilogram case.

To activate the SMART-PORT simply insert the power plug into an outlet. The integrated electric micromotor with LED light and the integrated electric scaler with LED light both have full performance capability.

Through the use of various contra-angles, speeds of between 0-200,000 rpm can be achieved. The high performance of the SMART-PORT is complemented by the highly efficient worldwide unique BPR Swiss suction (300 I/min!).

In addition to its unparalleled performance, the SMART-PORT offers all-around convenience. The clever design allows for quick and easy assembly and disassembly of the unit. In addition, the SMART-PORT has a remarkably quiet oil-free compressor and convenient externally removable instrument hoses.

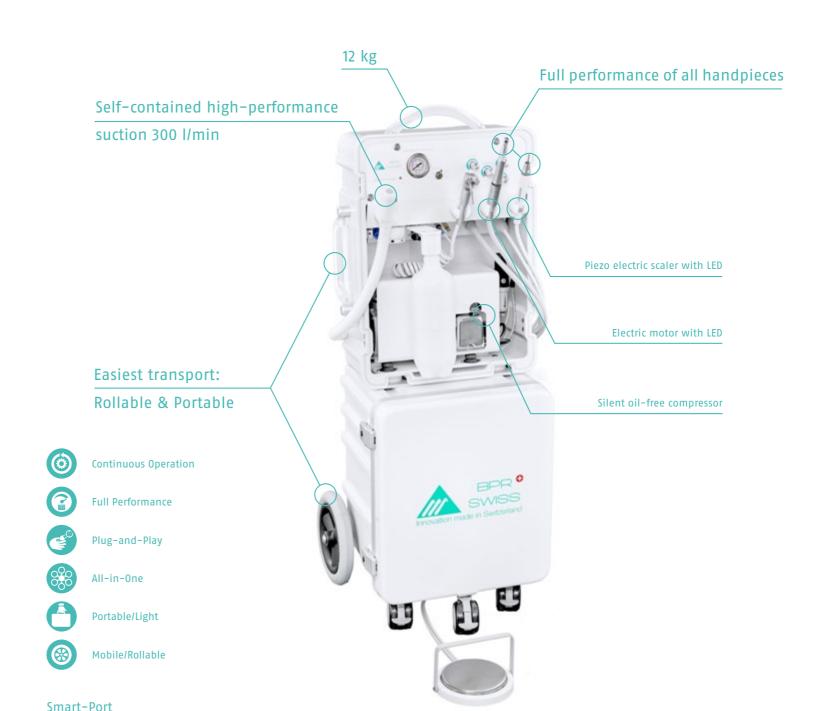



4 drawers for instruments/ consumables Separation red and green zone



On a single trolley transport of unit/instruments/consumables



Rollable and portable: In case of obstacles or for stowage the unit is undocked

Plug-and-play portable dental units

27